

00/23/20

2.24.00



JC525 U.S. PTO

Box: Patent Application  
 ASSISTANT COMMISSIONER FOR PATENTS  
 Washington, D.C. 20231

FORM PTO-1082  
 Case Docket No. PD-0329  
 Date: February 23, 2000  
 Express Mail Label No. EL 469 835 819 US

Dear Sir:

Transmitted herewith for filing is the patent application of  
 Inventor(s): **JOHN J. MASTROTOTARO; TODD M. GROSS; & JOHN J. SHIN**  
 For: **GLUCOSE MONITOR CALIBRATION METHODS**



Enclosed are:

- ☒ 39 page(s) application including cover sheet, 42 claims, and 1 page abstract  
☒ 10 Sheet(s) of informal drawing(s) (Figs. 1-15)  
☐ An assignment of the invention to MINIMED INC.  
☐ Declaration  
☐ Power of Attorney by Assignee

CALCULATION OF FEES						
ITEM		NO. OF CLAIMS FILED MINUS BASE*	NO. OF CLAIMS OVER BASE	X SM/LG ENTITY FEE	\$ AMOUNT	\$ FEE
A	TOTAL CLAIMS FEE	42- 20* =	22	X \$9 or X \$18	\$396	
B	INDEPENDENT CLAIMS FEE**	6- 3* =	3	X \$39 or X 78	\$234	
C	SUBTOTAL - ADDITIONAL CLAIMS FEE (ADD FINAL COLUMN IN LINES A + B)					\$630
D	SMALL ENTITY FEE = \$130 MULTIPLE-DEPENDENT CLAIMS FEE LARGE ENTITY FEE = \$260					\$
E	SMALL ENTITY FEE = \$380 BASIC FEE*LARGE ENTITY FEE = \$760					\$760
F	TOTAL FILING FEE (ADD TOTALS FOR LINES C, D, AND E)					\$1,390
F	ASSIGNMENT RECORDING FEE       \$ 40					\$40
	**LIST INDEPENDENT CLAIMS 1, 13, 34, 35, 38, 42					

- ☒ Please charge my Deposit Account No. 50-0621 the amount of \$1,430 **A copy of this letter is enclosed.**
- ☐ A check in the amount of \$ to cover the filing fee is enclosed.
- ☐ A check in the amount of \$ to cover Assignment Recordation fee is enclosed.

- ☒ The Commissioner is hereby authorized to charge payment of the following fees associated with this communication or credit any overpayment to Deposit Account No. 50-0621
- ☒ Any additional filing fees required under 37 CFR 1.16.
- ☒ Any patent application processing fees under 37 CFR 1.17.



02-24-00

A



Express Mail No. EL 469 835 819 US

PATENT  
PD-0329

**IN THE UNITED STATES PATENT AND TRADEMARK OFFICE**

In re Application of:	)	
John J. Mastrototaro et al.	)	Group Art Unit: unknown
Serial No.: unknown	)	
Filed: February 23, 2000	)	Examiner: unknown
For: GLUCOSE MONITOR	)	
CALIBRATION METHODS	)	
_____	)	

**CERTIFICATE OF U.S. EXPRESS MAIL**

Box: Patent Application  
Assistant Commissioner for Patents  
Washington, D.C. 20231

Sir:

Express Mail Mailing Label No. EL 469 835 819 US  
Date of Deposit: February 23, 2000

I hereby certify that a 49 page(s) patent application including cover sheet, 37 pages of specification, 42 claims, 1 page abstract, and 10 pages of informal drawings (Figs. 1 - 15), transmittal letter, and return postcard are being deposited using the "Express Mail Post Office to Addressee" service of the United States Postal Service under 37 C.F.R. § 1.10 on the above date and addressed to Assistant Commissioner for Patents, Washington, D.C. 20231.

February 23, 2000

date of deposit

Nestor Atienza

name of person mailing papers

Nestor Atienza

Signature

MiniMed Inc.  
12744 San Fernando Road  
Sylmar, CA 91342  
Tel.: (818) 362-2358 ext. 3313  
Fax.: (818) 367-1458

[illegible]

EL469835819US

TITLE

## Glucose Monitor Calibration Methods

RELATED APPLICATIONS

5           This application claims priority on U.S. Provisional Application Serial No. 60/121,584, filed February 25, 1999 and entitled "Glucose Monitor Calibration Techniques Using Linear Regression", which is herein incorporated by reference in its entirety.

FIELD OF THE INVENTION

10           This invention relates to glucose monitor systems and, in particular embodiments, to calibration methods for glucose monitoring systems.

BACKGROUND OF THE INVENTION

15           Over the years, body characteristics have been determined by obtaining a sample of bodily fluid. For example, diabetics often test for blood glucose levels. Traditional blood glucose determinations have utilized a painful finger prick using a lancet to withdraw a small blood sample. This results in discomfort from the lancet as it contacts nerves in the subcutaneous tissue. The pain of lancing and the cumulative discomfort from multiple needle pricks is a strong reason why patients fail to comply with a medical testing regimen used to  
20           determine a change in a body characteristic over a period of time. Although non-invasive systems have been proposed, or are in development, none to date have been commercialized that are effective and provide accurate results. In addition, all of these systems are designed to provide data at discrete points and do not provide continuous data to show the variations in the characteristic between testing times.

25           A variety of implantable electrochemical sensors have been developed for detecting and/or quantifying specific agents or compositions in a patient's blood. For instance, glucose

sensors are being developed for use in obtaining an indication of blood glucose levels in a diabetic patient. Such readings are useful in monitoring and/or adjusting a treatment regimen which typically includes the regular administration of insulin to the patient. Thus, blood glucose readings improve medical therapies with semi-automated medication infusion pumps of the external type, as generally described in U.S. Patent Nos. 4,562,751; 4,678,408; and 4,685,903; or automated implantable medication infusion pumps, as generally described in U.S. Patent No. 4,573,994, which are herein incorporated by reference. Typical thin film sensors are described in commonly assigned U.S. Patent Nos. 5,390,671; 5,391,250; 5,482,473; and 5,586,553 which are incorporated by reference herein. See also U.S. Patent No. 5,299,571.

#### SUMMARY OF THE DISCLOSURE

It is an object of an embodiment of the present invention to provide an improved glucose monitor system and method, which obviates for practical purposes, the above mentioned limitations.

According to an embodiment of the invention, a method of calibrating glucose monitor data includes obtaining glucose monitor data at predetermined intervals over a period of time. It also includes obtaining at least two reference glucose values from a reference source that correspond with the glucose monitor data obtained at the predetermined intervals. Additionally, calculating calibration characteristics using the at least two reference values and the corresponding glucose monitor data to regress the obtained glucose monitor data is included. And calibrating the obtained glucose monitor data using the calibration characteristics is included. In preferred embodiments, the reference source is a blood glucose meter, and the at least two reference glucose values are obtained from blood tests. In additional embodiments, the calculation of the calibration characteristics is obtained using linear regression, and in particular embodiments, using least squares linear regression. Alternatively, the calculation of the calibration characteristics is obtained using non-linear regression or a non-regression technique.

In particular embodiments, the predetermined period of time is a 24 hour period, and the predetermined intervals are 5 minute intervals. Further embodiments may include the step of shifting the data by a predetermined time factor, such as for example, ten minutes. Preferably, the calibration is performed while obtaining glucose monitor data. However, alternative  
5      embodiments may perform the calibration on glucose monitor data that has been collected for post processing by another processing device.

According to an embodiment of the invention, a method of calibrating glucose monitor data includes obtaining glucose monitor data at a predetermined memory storage rate. Also included is obtaining at least one blood glucose reference reading from a blood glucose  
10      measuring device that corresponds with at least one glucose monitor data point obtained at the predetermined memory storage rate. Calculating a calibration factor using the at least one blood glucose reference reading and the corresponding at least one glucose monitor data point is included. And calibrating the obtained glucose monitor data using the calibration factor is included. In preferred embodiments, after a first calibration factor is calculated, at least one  
15      previous calibration factor is used with at least one blood glucose reference reading from a blood glucose measuring device and its at least one corresponding glucose monitor data point to calculate a calibration factor. In additional embodiments, at least two blood glucose reference readings are used for calibration. In further embodiments, the calculation of the calibration factor is obtained using linear regression, and in particular least squares linear regression.  
20      Alternatively, calculation of the calibration factor uses non-linear regression or a non-regression technique

In particular embodiments, the calibration factor is applied to glucose monitor data obtained before a last blood glucose reference reading from a blood glucose measuring device that corresponds with at least one glucose monitor data point obtained at a predetermined  
25      memory storage rate is used to calculate the calibration factor. Alternatively, the calibration factor is applied to glucose monitor data obtained after the last blood glucose reference reading

from a blood glucose measuring device that is used to calculate the calibration factor.

In particular embodiments, the predetermined memory storage rate is once every 5 minutes. And the glucose monitor data that is obtained at a predetermined memory storage rate is the result of utilizing at least 2 sample values sampled from a glucose sensor at a rate faster than the memory storage rate.

In preferred embodiments, at least one blood glucose reference reading from a blood glucose measuring device is obtained during a predetermined calibration period, and a calibration factor is calculated using those readings after every predetermined calibration period. In particular embodiments, the predetermined calibration period is 24 hours. In further preferred embodiments, a predetermined time shift is used to temporally correlate the at least one blood glucose reference reading from a blood glucose measuring device with the at least one glucose monitor data point obtained at the predetermined memory storage rate. In particular embodiments, the predetermined time shift is ten minutes.

In particular embodiments, one or more calculations for calculating a first calibration factor is different from the one or more calculations for calculating subsequent calibration factors. In other particular embodiments, the calculation for calculating a first calibration factor uses a single-point calibration equation. In further particular embodiments, the single-point calibration equation includes an offset value. In other particular embodiments, the one or more calculations for calculating a calibration factor other than the first calibration factor uses a linear regression calibration equation, a non-linear regression calibration equation, or a non-regression technique.

According to an embodiment of the invention, a method of calibrating glucose monitor data includes obtaining glucose monitor data. It also includes obtaining from another blood glucose measuring device at least one blood glucose reference reading that is temporally associated with at least one glucose monitor data reading. Determining a calibration equation using the at least one blood glucose reference reading and the corresponding at least one glucose



monitor data reading is also included. And calibrating the glucose monitor data using the calibration equation is included.

According to another embodiment of the invention, a method of calibrating body characteristic monitor data includes obtaining body characteristic monitor data. It also includes obtaining from another characteristic measuring device at least one characteristic reference reading that is temporally associated with at least one characteristic monitor data point. Calculating calibration characteristics using the at least one characteristic reference reading and the corresponding at least one characteristic monitor data point is included. And calibrating the obtained characteristic monitor data using the calibration characteristics is included. In particular embodiments, at least two body characteristic reference readings are used for calculating the calibration characteristics. In particular embodiments, the calculation for calculating the calibration characteristics is a linear regression calculation.

According to additional embodiments of the invention, an apparatus for calibrating glucose monitor data includes a glucose monitor, glucose sensor, a blood glucose meter and a processor. The glucose monitor includes a glucose monitor memory for storing glucose monitor data. The glucose sensor is electronically coupled to the glucose monitor to supply the glucose monitor data. The blood glucose measuring device provides at least one blood glucose reference reading that is temporally associated with at least one glucose monitor data point. And the processor includes software to calculate calibration characteristics using the at least one blood glucose reference reading that is temporally associated with at least one glucose monitor data point, and the processor applies the calibration characteristics to the glucose monitor data. In particular embodiments, the at least one blood glucose reading is entered into the glucose monitor. In particular embodiments, the glucose monitor includes the processor, or alternatively, the processor is in a separate device that receives glucose monitor data from the glucose monitor.

In other embodiments of the invention, an apparatus for calibrating glucose monitor data includes means for obtaining glucose monitor data. It also includes means for obtaining from

another blood glucose measuring device at least one blood glucose reference reading that is temporally associated with at least one glucose monitor data reading. Means for calculating a calibration equation using the at least one blood glucose reference reading and the corresponding at least one glucose monitor data reading is included. And means for calibrating the glucose monitor data using the calibration equation is also included.

Other features and advantages of the invention will become apparent from the following detailed description, taken in conjunction with the accompanying drawings which illustrate, by way of example, various features of embodiments of the invention.

#### BRIEF DESCRIPTION OF THE DRAWINGS

A detailed description of embodiments of the invention will be made with reference to the accompanying drawings, wherein like numerals designate corresponding parts in the several figures.

Fig. 1 is a perspective view illustrating a subcutaneous glucose sensor insertion set and glucose monitor device in accordance with an embodiment of the present invention;

Fig. 2 is a cross-sectional view of the sensor set and glucose monitor device as shown along the line 2-2 of Fig. 1;

Fig. 3 is a cross-sectional view of a slotted insertion needle used in the insertion set of Figs. 1 and 2;

Fig. 4 is a cross-sectional view as shown along line 4-4 of Fig. 3;

Fig. 5 is a cross-sectional view as shown along line 5-5 of Fig. 3;

Fig. 6 is a partial cross-sectional view corresponding generally with the encircled region 6 of Fig. 2;

Fig. 7 is a cross-sectional view as shown along line 7-7 of Fig. 2.

Figs. 8(a-c) are diagrams showing a relationship between sampled values, interval values and memory storage values.

Fig. 9 is a chart showing clipping limits.

Fig. 10 is a sample computer screen image of a post processor analysis of glucose monitor data.

Fig. 11 is a chart illustrating the pairing of a blood glucose reference reading with glucose monitor data.

Fig. 12 is a chart illustrating an example of a single-point calibration.

Fig. 13 is a block diagram of a single-point calibration technique.

Fig. 14 is a chart illustrating an example of a linear regression calibration.

Fig. 15 is a block diagram of a linear regression calibration technique.

#### DETAILED DESCRIPTION OF THE PREFERRED EMBODIMENTS

As shown in the drawings for purposes of illustration, the invention is embodied in calibration methods for a glucose monitor that is coupled to a sensor set to provide continuous data recording of readings of glucose levels from a sensor for a period of time. In preferred embodiments of the present invention, the sensor and monitor are a glucose sensor and a glucose monitor for determining glucose levels in the blood and/or bodily fluids of a user. However, it will be recognized that further embodiments of the invention may be used to determine the levels of other body characteristics including analytes or agents, compounds or compositions, such as hormones, cholesterol, medications concentrations, viral loads (e.g., HIV), bacterial levels, or the like. The glucose sensor is primarily adapted for use in subcutaneous human tissue. However, in still further embodiments, one or more sensors may be placed in other tissue types, such as muscle, lymph, organ tissue, veins, arteries or the like, and used in animal tissue to measure body characteristics. Embodiments may record readings from the sensor on an intermittent, periodic, on-demand, continuous, or analog basis.

Figs. 1-7 illustrate a glucose monitor system 1 for use with the calibration methods. The glucose monitor system 1, in accordance with a preferred embodiment of the present invention,

includes a subcutaneous glucose sensor set 10 and a glucose monitor 100. In preferred embodiments, the glucose monitor 100 is of the type described in U.S. Patent Application Serial No. 60/121,664, filed on February 25, 1999, entitled "Glucose Monitor System", which is herein incorporated by reference. In alternative embodiments, the glucose monitor is of the type  
5 described in U.S. Patent Application Serial No. 09/377,472, filed August 19, 1999, entitled "Telemetered Characteristic Monitor System And Method Of Using The Same", which is incorporated by reference herein.

Preferably, the glucose monitor 100 is worn by the user and is connected to a surface mounted glucose sensor set 10 that is attached to a user's body by an electrically conductive  
10 cable 102, of the type described in U.S. Patent Application Serial No. 60/121,656, filed on February 25, 1999, entitled "Test Plug and Cable for a Glucose Monitor", which is incorporated by reference herein. In preferred embodiments, the sensor interface may be configured in the form of a jack to accept different types of cables that provide adaptability of the glucose monitor 100 to work with different types of subcutaneous glucose sensors and/or glucose sensors placed  
15 in different locations of the user's body. However, in alternative embodiments, the sensor interface is permanently connected to the cable 102. In additional alternative embodiments, a characteristic monitor is connected to one or more sensor sets to record data of one or more body characteristics from one or more locations on or in the user's body.

The glucose sensor set 10 is of the type described in U.S. Patent Application Serial No.  
20 60/121,655, filed on February 25, 1999, entitled "Glucose Sensor Set", or U.S. Patent Serial No. 08/871,831, filed on June 9, 1997, entitled "Insertion Set For A Transcutaneous Sensor", which are incorporated by reference herein. The glucose sensor 12, of the type described in U.S. Patent Application 29/101,218, filed on February 25, 1999, entitled "Glucose Sensor", or described in commonly assigned U.S. Patent Nos. 5,390,671; 5,391,250; 5,482,473; and 5,586,553 which are  
25 incorporated by reference herein; extends from the glucose sensor set 10 into the user's body with electrodes 20 of the glucose sensor 12 terminating in the user's subcutaneous tissue. See

also U.S. Patent No. 5,299,571. However, in alternative embodiments, the glucose sensor 12 may use other types of sensors, such as chemical based, optical based, or the like. In further alternative embodiments, the sensors may be of a type that is used on the external surface of the skin or placed below the skin layer of the user for detecting body characteristics.

5           The glucose monitor 100 generally includes the capability to record and store data as it is received from the glucose sensor 12, and includes either a data port (not shown) or wireless transmitter and/or receiver (also not shown) for transferring data to and/or from a data processor 200 such as a computer, communication station, a dedicated processor designed specifically to work with the glucose monitor, or the like. The glucose monitor is generally of the type  
10 described in U.S. Patent Application Serial No. 09/377,472, filed August 19, 1999, entitled “Telemetered Characteristic Monitor System And Method of Using The Same”, which is incorporated by reference herein.

          Preferably, the glucose monitor system 1 minimizes inconvenience by separating complicated monitoring process electronics into two separate devices; the glucose monitor 100,  
15 which attaches to the glucose sensor set 10; and the data processor 200, which contains the software and programming instructions to download and evaluate data recorded by the glucose monitor 100. In addition, the use of multiple components (e.g., glucose monitor 100 and data processor 200) facilitates upgrades or replacements, since one module, or the other, can be modified, re-programmed, or replaced without requiring complete replacement of the monitor  
20 system 1. Further, the use of multiple components can improve the economics of manufacturing, since some components may require replacement on a more frequent basis, sizing requirements may be different for each module, different assembly environment requirements, and modifications can be made without affecting the other components.

          The glucose monitor 100 takes raw glucose sensor data from the glucose sensor 12 and  
25 assesses it during real-time and/or stores it for later processing or downloading to the data processor 200, which in turn analyzes, displays, and logs the received data. The data processor

200 utilizes the recorded data from the glucose monitor 100 to analyze and review the blood glucose history. In particular embodiments, the glucose monitor 100 is placed into a com-station which facilitates downloading data to a personal computer for presentation to a physician. A software is used to download the data, create a data file, calibrate the data, and display the data in various formats including charts, forms, reports, graphs, tables, lists, and the like. In further  
5       embodiments, the glucose monitor system 1 may be used in a hospital environment or the like.

In alternative embodiments, the glucose monitor includes at least portions of the software described as contained within the data processor 200 above. The glucose monitor might contain the necessary software to calibrate glucose sensor signals, display a real-time blood glucose  
10       value, show blood glucose trends, activate alarms and the like. A glucose monitor with these added capabilities is useful for patients that might benefit from real-time observations of their blood glucose characteristics even while they're not in close proximity to a computer, communication device or dedicated independent data processor.

As shown in Fig. 2, the data processor 200, may include a display 214 that is used to  
15       display the calculated results of the raw glucose sensor data received via a download from the glucose monitor 100. The results and information displayed includes, but is not limited to, trending information of the characteristic (e.g., rate of change of glucose), graphs of historical data, average characteristic levels (e.g., glucose), stabilization and calibration information, raw data, tables (showing raw data correlated with the date, time, sample number, corresponding  
20       blood glucose level, alarm messages, and more), and the like. Alternative embodiments include the ability to scroll through the data. The display 214 may also be used with buttons (not shown) on the data processor 200, computer, communication station, characteristic monitor, or the like, to program or update data. In preferred embodiments, the glucose monitor 100 includes a display 132 to assist the user in programming the glucose monitor 100, entering data, stabilizing,  
25       calibrating, downloading data, or the like.

Still further embodiments of the present invention may include one or more buttons 122, 124, 126 and 128 on the glucose monitor 100 to program the monitor 100, to record data, insert flags to correlate data with external events for later analysis, input calibration values, or the like. In addition, the glucose monitor 100 may include an on/off button 130 for compliance with

5 safety standards and regulations to temporarily suspend transmissions or recording. The glucose monitor 100 may also be combined with other medical devices to accept other patient data through a common data network and/or telemetry system. The glucose monitor 100 may be combined with a blood glucose meter to directly import or correlate glucose calibration reference values such as described in U.S. Patent Application Serial No. 09/334,996, filed June 17, 1999,

10 entitled "Characteristic Monitor With A Characteristic Meter and Method Of Using The Same", which is incorporated by reference herein. The glucose monitor 100 may also be combined with semi-automated medication infusion pumps of the external type, as generally described in U.S. Patent Nos. 4,562,751; 4,678,408; and 4,685,903; or automated implantable medication infusion pumps, as generally described in U.S. Patent No. 4,573,994, which are herein incorporated by

15 reference. The glucose monitor 100 may record data from the infusion pumps and/or may process data from both the glucose sensor 12 and an infusion pump to establish a closed loop system to control the infusion pump based on glucose sensor measurements. In other embodiments, other body characteristics are monitored, and the monitor may be used to provide feedback in a closed loop system to control a drug delivery rate. In further alternative

20 embodiments, the glucose monitor 100 can be combined with a glucose sensor set 10 as a single unit.

Glucose sensors are replaced periodically to avoid infection, decaying enzyme coating and therefore sensor sensitivity, deoxidization of the electrodes, and the like. The user will disconnect the glucose sensor set 10 from the cable 102 and glucose monitor 100. A needle 14 is

25 used to install another glucose sensor set 10 and then the needle 14 is removed. Further description of the needle 14 and the sensor set 10 are found in U.S. Patent No. 5,586,553, entitled

“Transcutaneous Sensor Insertion Set”; co-pending U.S. Patent Application Serial No. 09/346,835, filed July 2, 1999, entitled “Insertion Set For A Transcutaneous Sensor”; and U.S. Patent No. 5,951,521, entitled “A Subcutaneous Implantable Sensor Set Having The Capability To Remove Or Deliver Fluids To An Insertion Site,” which are herein incorporated by reference.

5           The user connects the connection portion 24 of the glucose sensor set 10 through the cable 102 to the glucose monitor 100, so that the glucose sensor 12 can then be used over a prolonged period of time. An initial reading may be downloaded from the glucose sensor set 10 and the glucose monitor 100 to the data processor 200, to verify proper operation of the glucose sensor 10 and the glucose monitor 100. In preferred embodiments, the glucose sensor set 10  
10       provides data to the glucose monitor 100 for one to seven days before replacement. Glucose sensors 12 may last in the user’s body for longer or shorter periods of time depending on the quality of the installation, cleanliness, the durability of the enzyme coating, deoxidization of the sensor, user’s comfort, and the like.

          After installation into the body, the glucose sensor 12 is initialized to achieve a steady  
15       state of operation before starting a calibration process. Preferably, power supplied by three series silver oxide 357 battery cells 110 in the glucose monitor 100 is used to speed the initialization of the glucose sensor 12. Alternatively, other power supplies may be used such as, different battery chemistries including lithium, alkaline, or the like, and different numbers of batteries, solar cells, a DC converter plugged into an AC socket (provided with proper electrical isolation), or the like.

20       The use of an initialization process can reduce the time for glucose sensor 12 stabilization from several hours to an hour or less. The preferred initialization procedure uses a two step process. First, a high voltage (preferably between 1.0-1.1 volts - although other voltages may be used) is applied between electrodes 20 of the sensor 12 for one to two minutes (although  
25       different time periods may be used) to allow the sensor 12 to stabilize. Then, a lower voltage (preferably between 0.5-0.6 volts - although other voltages may be used) is applied for the remainder of the initialization process (typically 58 minutes or less). Other



stabilization/initialization procedures using differing currents, currents and voltages, different numbers of steps, or the like, may be used. Other embodiments may omit the initialization/stabilization process, if not required by the body characteristic sensor or if timing is not a factor. Alternatively, the characteristic monitor or the data processor 200 may apply an algorithm to the sensor data to determine when initial transients are sufficiently diminished and the sensor is at a significantly stable state to begin calibration.

In preferred embodiments, data is not considered valid until a sensor initialization event flag (ESI) is set in the data indicating that stabilization is complete. Preferably, stabilization is complete after 60 minutes or when a user enters a sensor initialization flag using one or more buttons on the glucose monitor 100. After stabilization/initialization is complete the glucose monitor 100 is calibrated to accurately interpret readings from the newly installed glucose sensor 12.

Beginning with the stabilization process, the glucose monitor 100 measures a continuous electrical current signal (ISIG) generated by the glucose sensor 12 relative to a concentration of glucose present in the subcutaneous tissue of the user's body. In preferred embodiments, the glucose monitor 100 samples the ISIG from the glucose sensor 12 at a sampling rate of once every 10 seconds, as shown in Fig. 8a-c. Examples of sampled values are labeled A - AD in Fig. 8a. At an interval rate of once per minute, the highest and lowest of the sampled values (shown in Fig. 8a as circled sampled values A, E, G, I, M, R, V, W, Y, and AB) are ignored, and the remaining 4 sampled values from an interval are averaged to create interval values (shown in Fig. 8b as values F', L', R', X', and AD'). At a glucose monitor memory storage rate of once every 5 minutes, the highest and lowest of the interval values (shown in Fig. 8b as values L' and X') are ignored and the remaining 3 interval values are averaged and stored in a glucose monitor memory as memory values (shown in Fig. 8c as point AD''). The memory values are retained in memory and may be downloaded to the data processor 200. The memory values are used to calibrate the glucose monitor 100 and/or the post processor 200 and to analyze blood glucose

levels. The sampling rate, interval rate and the memory storage rate may be varied as necessary to capture data with sufficient resolution to observe transients or other changes in the data depending on the rate at which sensor values can change, which is affected by the sensor sensitivity, the body characteristic being measured, the physical status of the user, and the like.

5 In other embodiments, all of the sampled values are included in the average calculations of memory storage values. In alternative embodiments, more or less sampled values or interval values are ignored depending on the signal noise, sensor stability, or other causes of undesired transient readings. Finally, in still other embodiments, all sampled values and/or interval values are stored in memory.

10 Clipping limits may be used to limit the signal magnitude variation from one value to the next thereby reducing the effects of extraneous data, outlying data points, or transients. In preferred embodiments, clipping limits are applied to the interval values. For instance, interval values that are above a maximum clipping limit or below a minimum clipping limit are replaced with the nearest clipping limit value.

15 In alternative embodiments, interval values that are outside of the clipping limits are ignored and not used to calculate the next memory storage value. In particular embodiments, the detection of interval values outside of the clipping limits is considered a calibration cancellation event. In further particular embodiments, more than one value must be deemed outside of clipping limits to constitute a calibration cancellation event. (Calibration cancellation events are  
20 discussed below).

In preferred embodiments, the clipping limits are shifted after each data point. The level that the clipping limits are set to is dependent on an acceptable amount of change from the previous interval value to the present interval value, which is affected by the sensor sensitivity, signal noise, signal drift, and the like. In preferred embodiments, the clipping limits are  
25 calculated based on the magnitude of the previous interval value. For example, for a previous interval value from 0 up to but not including 15 Nano-Amps, the clipping limits are set at plus

and minus 0.5 Nano-Amps about the previous interval value. For a previous interval value from 15 up to but not including 25 Nano-Amps, the clipping limits are set at plus and minus 3% of the previous interval value, about the previous interval value. For a previous interval value from 25 up to but not including 50 Nano-Amps, the clipping limits are set at plus and minus 2% of the previous interval value, about the previous interval value. And for a previous interval value of 50 Nano-Amps and greater, the clipping limits are set at plus and minus 1% about the previous interval value. In alternative embodiments, different clipping limits may be used.

Fig. 9 shows a typical clipping limit example in which a previous interval value 500, associated with interval N-1, has a magnitude of 13.0 Nano-Amps, which is less than 15.0 Nano-Amps. Therefore, the maximum clipping limit 502 for the present interval value 506 is set at 13.5 Nano-Amps, which is 0.5 Nano-Amps greater than the magnitude of the previous interval value 500. And the minimum clipping limit 504 is set at 12.5 Nano-Amps which is 0.5 Nano-Amps below the previous interval value 500. The present interval value 506, associated with interval N, is between the maximum clipping limit 502 and the minimum clipping limit 504 and is therefore acceptable.

In another example shown in Fig. 9, the present interval value 508, associated with interval M, has a value of 25.0 Nano-Amps which is outside of the clipping limit 514 and will therefore be clipped. The previous interval value 510, associated with interval M-1, is 26.0 Nano-Amps, which is included in the range from 25.0 up to but not including 50.0 Nano-Amps as discussed above. Therefore the clipping limits are  $\pm 2\%$ . The maximum clipping limit 512 is 2% greater than the previous interval value 510,

$$26.0 + 26.0 * 0.02 = 26.5 \text{ Nano-Amps.}$$

Similarly the minimum clipping limit 514 is 2% less than the previous interval value 510,

$$26.0 - 26.0 * 0.02 = 25.5 \text{ Nano-Amps.}$$

Since the present interval value 508 of 25.0 Nano-Amps is less than the minimum clipping limit 514 of 25.5 Nano-Amps, it will be clipped, and 25.5 Nano-Amps will be used in place of 25.0

Nano-Amps to calculate a memory storage value. For further illustration, Fig. 8 shows interval value  $R'$ , which is calculated by averaging sampled values  $N$  through  $Q$ , is outside of the clipping limits 412 and 414, which result from the previous interval value  $L'$ . Therefore, the magnitude of interval value  $R'$  is not used to calculate memory value  $AD''$ , instead  $R''$ , which is the magnitude of the minimum clipping limit 414, is used.

In other embodiments, the clipping limits may be a smaller or larger number of Nano-Amps or a smaller or larger percentage of the previous interval value based on the sensor characteristics mentioned above. Alternatively, the clipping limits are calculated as plus or minus the same percent change from every previous interval value. Other algorithms use several interval values to extrapolate the next interval value and set the clipping limits to a percentage higher and lower than the next anticipated interval value. In further alternatives, clipping may be applied to the sampled values, interval values, memory values, calculated glucose values, estimated values of a measured characteristic, or any combination of the values.

In preferred embodiments, all interval values are compared to an out-of-range limit of 200 Nano-Amps. If three consecutive interval values are equal to or exceed the out-of-range limit, the sensor sensitivity is deemed to be too high and an alarm is activated to notify the user that re-calibration is required or the sensor may need replacing. In alternative embodiments, the out-of-range limit is set at higher or lower values depending on the range of sensor sensitivities, the expected working life of the sensor, the range of acceptable measurements, and the like. In particular embodiments, the out-of range limit is applied to the sampled values. In other embodiments, the out-of-range limit is applied to the memory storage values.

In preferred embodiments, unstable signal alarm limits are set to detect when memory storage values change too much from one to another. The signal alarm limits are established similarly to the clipping limits described above for the interval values, but allow for a larger change in value since there is more time between memory storage values than between interval values. Re-calibration or replacement of the glucose sensor 12 is required once an unstable

signal alarm is activated. In essence, the glucose monitor 100 has detected too much noise in the ISIG from the glucose sensor 12.

Each memory storage value is considered valid (Valid ISIG value) unless one of the following calibration cancellation events occurs: an unstable signal alarm (as discussed above), a sensor initialization event (as discussed above), a sensor disconnect alarm, a power on/off event, an out-of-range alarm (as discussed above), or a calibration error alarm. Only Valid ISIG values are used to calculate blood glucose levels by the glucose monitor 100 or post processor 200, as shown in Fig. 10. Once a calibration cancellation event occurs, the successive memory storage values are not valid, and therefore are not used to calculate blood glucose, until the glucose monitor 100 or post processor 200 is re-calibrated. Fig. 10 shows an explanatory computer screen in which cell P3 indicates a sensor disconnect alarm with the abbreviation "SeDi". As shown, blood glucose values do not appear in column K, titled "Sensor Value", and Valid ISIG values do not appear in column J until after the sensor is initialized, as indicated by the "ESI" flag in cell N17. One exception however, is the power on/off event. If the glucose monitor 100 is turned off for a short enough period of time, generally up to 30 minutes, the memory storage values are considered Valid ISIG values as soon as the power is turned back on. If the power is off for longer than 30 minutes, the glucose monitor must be re-calibrated before ISIG values are considered valid. Alternatively, the power may be off 30 minutes up to indefinitely and once the power is restored, the memory storage values are Valid ISIG values. The sensor disconnect alarm is activated when the glucose monitor 100 does not detect a signal. In preferred embodiments, when 2 or more out of 5 interval values collected within a given memory storage rate are less than 1.0 Nano-Amp, the disconnect alarm is triggered. In alternative embodiments, more or less values need be below a particular amperage to trigger the disconnect alarm depending of the acceptable range or sensor readings and the stability of the sensor signal. The remaining two calibration cancellation events, the calibration error and an alternative embodiment for the out-of-range alarm, are discussed in conjunction with the calibration process

below.

Preferred embodiments are directed to calibration techniques that are used by either glucose monitors 100 during real-time measurements of one or more signals from the glucose sensor 12, or post processors 200 during post-processing of data that has been previously  
 5 recorded and downloaded (as shown in Fig. 10).

To calibrate the glucose monitor 100, the calibration factor called a sensitivity ratio (SR) (blood glucose level/Valid ISIG value) is calculated for a particular glucose sensor 12. The SR is a calibration factor used to convert the Valid ISIG value (Nano-Amps) into a blood glucose level (mg/dℓ or mmol/ℓ). In alternative embodiments, the units for the SR may vary depending  
 10 on the type of signal available from the sensor (frequency, amplitude, phase shift, delta, current, voltage, impedance, capacitance, flux, and the like), the magnitude of the signals, the units to express the characteristic being monitored, or the like.

In preferred embodiments, the user obtains a blood glucose reference reading from a common glucose meter, or another blood glucose measuring device, and immediately enters the  
 15 blood glucose reference reading into the glucose monitor 100. The blood glucose reference reading is assumed to be accurate and is used as a reference for calibration. The glucose monitor 100, or a post processor 200, must temporally correlate the blood glucose reference reading with a Valid ISIG value to establish a paired calibration data point. Since the glucose level in the interstitial body fluid tends to lag behind the blood glucose level, the glucose monitor 100 or post  
 20 processor 200 applies a delay time and then pairs the blood glucose reference reading with a Valid ISIG value as shown in Fig. 11. In preferred embodiments, an empirically derived 10 minute delay is used. Since Valid ISIG values are averaged and stored every 5 minutes, the glucose monitor 100 correlates the blood glucose reference reading with the third Valid ISIG stored in memory after the blood glucose reference reading is entered (resulting in an effective  
 25 delay of 10 to 15 minutes). Fig. 11 illustrates an example, in which a blood glucose reference reading 600 of 90 mg/dℓ is entered into the glucose monitor 100 at 127 minutes. The next Valid

ISIG value 602 is stored at 130 minutes. Given a 10 minute delay, the glucose reference reading 600 is paired with the Valid ISIG value 604 which is stored at 140 minutes with a value of 30 Nano-amps. Note that two numbers are needed to establish one paired calibration data point, a blood glucose reference reading and a Valid ISIG.

5           Other delay times may be used depending on the user's metabolism, the response time of the sensor, the delay time required for the glucose meter to calculate a reading and for the reading to be entered into the glucose monitor 100, the type of analyte being measured, the tissue that the sensor is placed into, environmental factors, whether the previous glucose Valid ISIG value (or the trend of the Valid ISIG values) was higher or lower than current Valid ISIG value, or the like. Once paired calibration data is available, the appropriate calibration process may be  
10           applied dependent on how many paired calibration data points are available since the last calibration, the total period of time that the glucose sensor 12 has been in use, and the number of times the glucose sensor 12 has been calibrated.

15           In preferred embodiments, blood glucose reference readings are entered into the glucose monitor 100 periodically through out each day of use. Preferably calibration is conducted immediately after the initialization/stabilization of a glucose sensor 12 and once a day thereafter. However, calibration may be conducted more or less often depending on whether a glucose sensor 12 has been replaced, whether a calibration cancellation event has occurred, the stability of the glucose sensor 12 sensitivity over time, or the like.

In preferred embodiments, blood glucose reference readings are collected several times per day but a new calibration factor is calculated only once per day. Therefore, typically more than one paired calibration data point is collected between calibrations. In alternative embodiments, the glucose monitor is calibrated every time a new paired calibration data point is collected.

Preferred embodiments use a single-point calibration technique (shown in a block diagram of Fig. 13) to calculate the SR when only one paired calibration data point is available, such as immediately after initialization/stabilization. And a modified linear regression technique (shown in a block diagram in Fig. 15) is used when two or more paired calibration data points are available. Particular embodiments use a single-point calibration technique whether or not more than one paired calibration data point is available.

A single-point calibration equation is based on the assumption that the Valid ISIG will be 0 when the blood glucose is 0. As shown in Fig. 12, a single paired calibration point 700 is used with the point (0,0) to establish a line 702. The slope of the line from the origin (0,0) and passing through the single paired calibration point 700 is the single-point sensitivity ratio (SPSR). The single-point calibration equation to calculate the calibration factor SPSR is as follows:

$$\text{SPSR} = \frac{\text{Blood Glucose Reference Reading}}{\text{Valid ISIG}}$$

Where SPSR = a Single-Point Sensitivity Ratio.

Therefore, the calibrated blood glucose level is,

$$\text{Blood Glucose Level} = \text{Valid ISIG} * \text{SPSR}.$$



As an example, using the values of 20.1 Nano-Amps and 102 mg/dℓ from the paired calibration data point shown in Fig. 12, the calculation of SPSR is:

$$\text{SPSR} = 102 / 20.1 = 5.07 \text{ mg/dℓ per Nano-amp.}$$

To continue the example, once the calibration is complete, given a glucose sensor reading of 15.0 Nano-Amps, the calculated blood glucose level is:

$$\text{Blood Glucose Level} = 15.0 * 5.07 = 76.1 \text{ mg/dℓ.}$$

Additionally, particular embodiments use an offset value in a calibration equation to compensate for the observation that more sensitive glucose sensors 12 (i.e. glucose sensors 12 that generate higher ISIG values compared to other glucose sensors 12 at the same blood glucose level, which result in lower SR values) often have a less linear performance at very high blood glucose levels when compared to glucose sensors 12 with lower sensitivity (and therefore relatively higher SR values). If the SPSR for a particular glucose sensor 12, as calculated above, is less than a sensitivity threshold value, then a modified SPSR (MSPSR) is calculated using an offset value included in a modified single-point calibration equation. In preferred embodiments, the threshold value is 7. When the initial calculation of the SPSR (shown above) is less than 7, an offset value of 3 is used to calculate the MSPSR. If the initial calculation of SPSR yields a value of 7 or greater, then the offset value is 0. Thus, the calibration factor (MSPSR) is calculated using the offset value in the modified single-point calibration equation, as follows:

$$\text{MSPSR} = \frac{\text{Blood Glucose Reference Reading}}{(\text{Valid ISIG} - \text{offset})}$$

Therefore, the calibrated blood glucose level is,

$$\text{Blood Glucose Level} = (\text{Valid ISIG} - \text{offset}) * \text{MSPSR}.$$

- 5 Continuing the above example since the SPSR is 5.07, which is less than 7, the sensitivity ratio is recalculated using the MSPSR equation as:

$$\text{MSPSR} = 102 / (20.1 - 3) = 5.96 \text{ mg/dl per Nano-amp.}$$

- 10 And given a glucose sensor reading of 15.0 Nano-Amps after calibration the calculated blood glucose level is:

$$\text{Blood Glucose Level} = (15.0 - 3) * 5.96 = 71.5 \text{ mg/dl}.$$

- 15 In another example, given a blood glucose reference reading of 95 from a typical blood glucose meter and a Valid ISIG value of 22.1, the resulting SPSR is  $95/22.1 = 4.3$ . Since  $SR < 7$ , the offset = 3. Therefore, the MSPSR is  $95 / [22.1 - 3] = 5.0$ . Note that when the SPSR is greater than or equal to 7 the offset value is 0 and therefore the  $\text{MSPSR} = \text{SPSR}$ .

- 20 In alternative embodiments, the offset value is eliminated from the equation for calculating the blood glucose value as follows:

$$\text{Blood Glucose Level} = \text{Valid ISIG} * \text{MSPSR}.$$

- 25 The threshold value of 7 and the associated offset of 3 have been empirically selected based on the characteristics observed from testing a particular type of glucose sensors 12, such as those described in U.S. Patent 5,391,250 entitled "Method of Fabricating Thin Film Sensors", and U.S.

Patent Application Serial Number 60/121,655 filed on February 25, 1999, entitled "Glucose Sensor Set", incorporated by reference herein. Other threshold values may be used in conjunction with other offset values to optimize the accuracy of the calculated MSPSR for various types of glucose sensors 12 and sensors used to detect other body characteristics. In fact, many threshold values may be used to select between many offset values. An example using two different threshold values (4 and 7) to select between three different offset values (5, 3 and 0) follows:

If the SPSR < 4, then use an offset value of 5, else  
 if  $4 \leq \text{SPSR} < 7$ , then use an offset value of 3, else  
 if  $\text{SPSR} \geq 7$  use an offset value of 0.

In preferred embodiments the MSPSR is compared to a valid sensitivity range to determine if the newly calculated MSPSR is reasonable. In order to identify potential system problems, a valid MSPSR range of 1.5 to 15 is employed. This range has been determined based upon valid glucose sensor sensitivity measurements made *in-vitro*. MSPSR values outside this range result in a calibration error alarm (CAL ERROR) to notify the user of a potential problem. Other valid sensitivity ranges may be applied depending on the types of sensors to be calibrated, the range of acceptable sensitivity levels for the various sensor types, the manufacturing consistency expected for the sensors, environmental conditions, how long the sensor has been in use, or the like.

Preferred embodiments augment the single-point calibration technique using a modified linear regression technique (shown in a block diagram in Fig. 15) when more than one paired calibration data point is available. As shown in Fig. 14, the paired calibration data points 800 are linearly regressed by a least squares method to calculate the best fit straight line 802 correlated with paired calibration data points 800. The slope of the line resulting from the linear regression is the linear regression sensitivity ratio (LRSR) used as the calibration factor to calibrate the glucose monitor 100. The linear regression calibration equation is as follows:

$$LRSR = \frac{\sum_{i=1}^N [X_i Y_i]}{\sum_{i=1}^N [X_i^2]}$$

Where  $X_i$  is the  $i$ th Valid ISIG value of paired calibration data points,

$Y_i$  is the  $i$ th Blood Glucose Reference Reading of paired calibration data points and,

5  $N$  is the total number of paired calibration data points used for calibration.

$i$  is the identification number of a particular paired calibration data point.

Therefore, the calibrated blood glucose level is,

$$10 \quad \text{Blood Glucose Level} = \text{Valid ISIG} * \text{LRSR}.$$

Note that this linear regression uses a fixed intercept of zero (in other words, when the Valid ISIG is 0 the blood glucose value is 0) and therefore the linear regression method estimates only one regression parameter, the slope. In alternative embodiments, other linear regression methods  
15 may be used that estimate additional regression parameters such as an offset value.

Additionally, particular embodiments use an offset value in a modified linear regression calibration equation. The purpose of the offset value, as described above for the single-point calibration, is to compensate for the observation that more sensitive glucose sensors 12 often have a less linear performance at very high blood glucose levels. If the LRSR for a particular glucose  
20 sensor 12, as calculated in the linear regression calibration equation above, is less than a sensitivity threshold value, then a modified linear regression sensitivity ratio (MLRSR) is calculated using an offset value included in a modified linear regression calibration equation. In preferred embodiments, the threshold value is 7. When the initial calculation of the LRSR is less than 7, an offset value of 3

is used to calculate the MLRSR. If the initial calculation of LRSR yields a value of 7 or greater, then the offset value is 0. Thus, the MLRSR is calculated using the offset value in the modified linear regression calibration equation, thus:

$$MLRSR = \frac{\sum_{i=1}^N [(X_i - offset)Y_i]}{\sum_{i=1}^N (X_i - offset)^2}$$

Therefore, the calibrated blood glucose level is,

$$\text{Blood Glucose Level} = (\text{Valid ISIG-offset}) * \text{MLRSR}.$$

Just as in the case of single-point calibration techniques described above, other threshold values may be used in conjunction with other offset values in the modified linear regression calibration equation to optimize the accuracy of the calculated MLRSR for various types of glucose sensors 12 and other characteristic sensors.

In preferred embodiments the MLRSR is compared to a valid sensitivity range to determine if the newly calculated MLRSR is reasonable. In order to identify potential system problems, a valid MLRSR range of 2.0 to 10.0 is employed. MLRSR values outside this range result in a calibration error alarm (CAL ERROR) to notify the user of a potential problem. As described above for the single-point calibration techniques, other valid sensitivity ranges may be applied.

In preferred embodiments, the glucose monitor data is linearly regressed over a 24 hour period (or window), and new sensitivity ratios are used for each 24 hour time period. In alternative embodiments, the time period may be reduced to only a few hours or enlarged to

cover the entire monitoring period with the glucose sensor (i.e., several days - or even weeks with implanted sensors). In further embodiments, the time window may be fixed at a predetermined size, such as 24 hours, 12 hours, 6 hours, or the like, and the window is moved along over the operational life of the sensor.

5 In particular embodiments, paired calibration data points from measurements taken before the last calibration may be used to calculate a new sensitivity ratio. For example, to calibrate the glucose monitor every 6 hours, a paired calibration data point is established every 6 hours. And the linear regression technique described above is executed using 4 paired calibration data points, the most recently acquired point and points from 6, 12 and 18 hours before.

10 Alternatively, the number of paired calibration data points used in the calibration may be as few as one or as large as the total number of paired calibration data points collected since the glucose sensor was installed. In alternative embodiments, the number of paired calibration data points used in a calibration equation may grow or shrink during the life of the glucose sensor due to glucose sensor anomalies.

15 In still other embodiments, the decay characteristics of the glucose sensor 12 over time may be factored into the equation to account for typical degradation characteristics of the glucose sensor 12 due to site characteristics, enzyme depletion, body movement, or the like. Considering these additional parameters in the calibration equation will more accurately tailor the calibration equation used by the glucose monitor 100 or post processor 200. In particular embodiments, 20 other parameters may be measured along with the blood glucose such as, temperature, pH, salinity, and the like. And these other parameters are used to calibrate the glucose sensor using non-linear techniques.

In alternative embodiments, the calibration factor may be calculated by first using a single-point technique to calculate the MSPSR for each paired calibration data point and then 25 averaging them together, either unweighted or weighted by temporal order or by elapsed time. In other alternative embodiments, other straight line curve fitting techniques may be used to

calculate a slope to be used as the SR. In additional alternative embodiments, other non-regressive curve fitting techniques may be used to generate equations that express the blood glucose level relative to the Valid ISIG. The equations may be polynomial, parabolic, hyperbolic, asymptotic, logarithmic, exponential, or the like. In these embodiments, the SR is not a single value (such as a slope) but rather an equation representing a curve that is used to convert the Valid ISIG from the glucose sensor 12 to a blood glucose value in the glucose monitor 100 or a post processor 200. This method is especially effective when using sensors that produce signals that are not linearly proportional to the characteristic being measured across the entire range of measurement.

As discussed, preferred embodiments utilize a least squares linear regression equation to calibrate the glucose monitor 100 or post-processor 200 to analyze the sensor data. However, alternative embodiments may utilize a multiple component linear regression, or equations with more variables than just the paired calibration data points, to account for additional calibration effecting parameters, such as environment, the individual user's characteristics, sensor lifetime, manufacturing characteristics (such as lot characteristics), deoxidization, enzyme concentration fluctuation or degradation, power supply variations, or the like. Still other alternative embodiments may utilize singular and multiple, non-linear regression techniques.

In preferred embodiments, after the first calibration is performed on a particular glucose sensor 12, subsequent calibrations employ a weighted average using a sensitivity ratio (SPSR, MSPSR, LRSR, or MLRSR) calculated from data collected since the last calibration, and previous sensitivity ratios calculated for previous calibrations. So the initial sensitivity ratio (SR1) is calculated immediately after initialization/stabilization using a paired calibration data point and is used by the glucose monitor 100 or the post processor 200 until the second sensitivity ratio (SR2) is calculated. The second sensitivity ratio (SR2) is an average of SR1 and the sensitivity ratio as calculated using the paired calibration data points since the initial calibration (SRday1). The equation is as follows:

$$SR2 = \frac{(SR1 + SRday1)}{2}$$

The third sensitivity ratio (SR3) is an average of SR2 and the sensitivity ratio as calculated using the paired calibration data points since the second calibration (SRday2). The equation is as follows:

$$SR3 = \frac{(SR2 + SRday2)}{2}$$

The sensitivity ratios for successive days use the same format, which is expressed below in generic terms:

$$SR_n = \frac{(SR_{(n-1)} + SRday_{(n-1)})}{2}$$

Where  $SR_n$  is the new sensitivity ratio calculated at the beginning of time period,  $n$ , using data from time period  $(n-1)$ , to be used by a real time glucose monitor 100, to convert Valid ISIGs to blood glucose readings throughout time period,  $n$ .

$SR_{(n-1)}$  is the previous sensitivity ratio calculated at the beginning of time period,  $n-1$ , using data from time period  $(n-2)$ .

$SRday_{(n-1)}$  is the sensitivity ratio calculated using paired calibration data points collected since the last calibration.

Alternatively, the previous sensitivity ratios may be ignored and the SR is calculated using only the paired calibration data points since the last calibration. Another alternative is to equally average all previous SRs with the latest SR calculated using only the paired calibration data points since the last calibration. In alternative embodiments, the paired calibration data points are used to establish an equation for a curve representing SR over time. The curve is then used to extrapolate SR to be used until the next paired calibration data point is entered.

In embodiments that use a post processor 200 to evaluate the sensitivity ratio, the



sensitivity ratio is calculated using paired calibration data points over a period of time since the last calibration and is not averaged with previous SRs. The sensitivity ratio for a period of time can then be applied to the same period of time over which the paired calibration data points were collected. This is more accurate than the real-time case described above for the glucose monitor 100 because, in the real-time case, sensitivity ratios from a previous time period must be used to calculate the blood glucose level in the present time period. If the sensitivity ratio has changed over time, the calculation of blood glucose using an old sensitivity ratio introduces an error.

In particular embodiments, once calibration is complete, Valid ISIG values are converted to blood glucose readings based on a particular version of the sensitivity ratio, and the resulting blood glucose readings are compared to an out-of-range limit. If the resulting calculated blood glucose level is greater than a maximum out-of-range limit of 200 mg/dℓ (or equivalently 3600 mmol/ℓ), the out-of-range alarm is activated. This is a calibration cancellation event, therefore, ISIG values are no longer valid once this alarm is activated. The blood glucose readings are either not calculated, or at least not considered reliable, until the glucose monitor 100 or post processor 200 is re-calibrated. The user is notified of the alarm and that re-calibration is needed.

In alternative embodiments, higher or lower maximum out-of-range limits may be used depending on the sensor characteristics, the characteristic being measured, the user's body characteristics, and the like. In particular embodiments, a minimum out-of-range limit may be used or both a maximum and a minimum out-of-range limits may be used. In other particular embodiments, the out-of-range limits do not cause the blood glucose readings to become invalid and/or re-calibration is not required; however, an alarm could still be provided. In additional particular embodiments, more than one ISIG value must exceed an out-of-range limit before an alarm is activated of a calibration cancellation event is triggered. The ISIG values that are out-of-range are not used to display a blood glucose value.

In alternative embodiments, calibration is conducted by injecting a fluid containing a known value of glucose into the site around the glucose sensor set 10, and then one or more

glucose sensor readings are sent to the glucose monitor 100. The readings are processed (filtered, smoothed, clipped, averaged, and the like) and used along with the known glucose value to calculate the SR for the glucose sensor 12. Particular alternative embodiments, use a glucose sensor set of the type described in U.S. Patent 5,951,521 entitled "A Subcutaneous Implantable Sensor Set Having the Capability To Remove Or Deliver Fluids To An Insertion Site".

In other alternative embodiments, the glucose sensor 12 is supplied with a vessel containing a solution with a known glucose concentration to be used as a reference, and the glucose sensor 12 is immersed into the reference glucose solution during calibration. The glucose sensor 12 may be shipped in the reference glucose solution. As described above, the glucose sensor readings are used to calculate a sensitivity ratio given the known glucose concentration of the solution.

In another alternative embodiment, the glucose sensors 12 are calibrated during the manufacturing process. Sensors from the same manufacturing lot, that have similar properties, are calibrated using a sampling of glucose sensors 12 from the population and a solution with a known glucose concentration. The sensitivity ratio is provided with the glucose sensor 12 and is entered into the glucose monitor 100 or the post processor 200 by the user or another individual.

While the description above refers to particular embodiments of the present invention, it will be understood that many modifications may be made without departing from the spirit thereof. The accompanying claims are intended to cover such modifications as would fall within the true scope and spirit of the present invention.

The presently disclosed embodiments are therefore to be considered in all respects as illustrative and not restrictive, the scope of the invention being indicated by the appended claims, rather than the foregoing description, and all changes which come within the meaning and range of equivalency of the claims are therefore intended to be embraced therein.

WHAT IS CLAIMED IS:

1. A method of calibrating glucose monitor data, the method comprising the steps of:  
obtaining glucose monitor data at predetermined intervals over a period of time;  
obtaining at least two reference glucose values from a reference source that correspond  
5 with the glucose monitor data obtained at the predetermined intervals;  
calculating calibration characteristics using the at least two reference values and the  
corresponding glucose monitor data to regress the obtained glucose monitor data; and  
calibrating the obtained glucose monitor data using the calibration characteristics.
- 10 2. The method of claim 1, wherein the reference source is a blood glucose meter, and  
wherein the at least two reference glucose values are obtained from blood tests.
3. The method of claim 1, wherein the calculation of the calibration characteristics is  
obtained using linear regression.
- 15 4. The method of claim 3, wherein the linear regression is least squares linear regression.
5. The method of claim 1, wherein the calculation of the calibration characteristics is  
obtained using non-linear regression.
- 20 6. The method of claim 1, wherein the calculation of the calibration characteristics is  
obtained using a non-regression technique.
7. The method of claim 1, wherein the predetermined period of time is a 24 hour period.
- 25 8. The method of claim 1, wherein the predetermined intervals are 5 minute intervals.

9. The method of claim 1, wherein the method further comprises the step of shifting the data by a predetermined time factor.

5 10. The method of claim 9, wherein the predetermined time factor is ten minutes.

11. The method of claim 1, wherein the calibration is performed while obtaining the glucose monitor data.

10 12. The method of claim 1, wherein the calibration is performed on glucose monitor data that has been collected for post processing by another processing device.

13. A method of calibrating glucose monitor data, the method comprising the steps of:  
obtaining glucose monitor data at a predetermined memory storage rate;

15 obtaining at least one blood glucose reference reading from a blood glucose measuring device that corresponds with at least one glucose monitor data point obtained at the predetermined memory storage rate;

calculating a calibration factor using the at least one blood glucose reference reading and the corresponding at least one glucose monitor data point; and

20 calibrating the obtained glucose monitor data using the calibration factor.

14. The method of claim 13, wherein after a first calibration factor is calculated, at least one previous calibration factor is used with the at least one blood glucose reference reading from a blood glucose measuring device that corresponds with at least one glucose monitor data point  
25 obtained since the at least one previous calibration factor was calculated, to calculate a new calibration factor.

15. The method of claim 13, wherein the at least one blood glucose reference reading is at least two blood glucose reference readings.

5 16. The method of claim 15, wherein a calculation of the calibration factor is obtained using linear regression.

17. The method of claim 16, wherein the linear regression is least squares linear regression.

10 18. The method of claim 13, wherein a calculation of the calibration factor is obtained using non-linear regression.

19. The method of claim 13, wherein a calculation of the calibration factor is obtained using a non-regression technique.

15 20. The method of claim 15, wherein the calibration factor is applied to glucose monitor data obtained before the last of the at least one blood glucose reference reading from a blood glucose measuring device that corresponds with at least one glucose monitor data point obtained at the predetermined memory storage rate used to calculate the calibration factor.

20 21. The method of claim 13, wherein the calibration factor is applied to glucose monitor data obtained after a last at least one blood glucose reference reading from a blood glucose measuring device that corresponds with at least one glucose monitor data point obtained at the predetermined memory storage rate used to calculate the calibration factor.

25

22. The method of claim 13, wherein the predetermined memory storage rate is once every 5 minutes.

23. The method of claim 13, wherein the glucose monitor data that is obtained at a  
5 predetermined memory storage rate is the result of utilizing at least 2 sample values sampled from a glucose sensor at a rate faster than the memory storage rate.

24. The method of claim 13, wherein at least one blood glucose reference reading from a blood glucose measuring device is obtained during a predetermined calibration period, and  
10 wherein a calibration factor is calculated using at least one blood glucose reference reading from a blood glucose measuring device after every predetermined calibration period.

25. The method of claim 24, wherein predetermined calibration period is a 24 hours.

26. The method of claim 13, wherein one or more calculations for calculating a first  
15 calibration factor is different than one or more calculations for calculating subsequent calibration factors.

27. The method of claim 26, wherein the calculation for calculating a first calibration factor  
20 uses a single-point calibration equation.

28. The method of claim 27, wherein the single-point calibration equation includes an offset value.

29. The method of claim 27, wherein the calculation for calculating the subsequent  
25 calibration factors uses a linear regression calibration equation.

30. The method of claim 27, wherein the calculation for calculating the subsequent calibration factors uses a non-linear regression calibration equation.

5 31. The method of claim 27, wherein the calculation for calculating the subsequent calibration factors uses a non-regression calibration equation.

32. The method of claim 13, wherein a predetermined time shift is used to temporally correlate the at least one blood glucose reference reading from a blood glucose measuring device  
10 with the at least one glucose monitor data point obtained at the predetermined memory storage rate.

33. The method of claim 13, wherein the predetermined time shift is ten minutes.

15 34. A method of calibrating glucose monitor data, the method comprising the steps of:  
obtaining glucose monitor data;  
obtaining from another blood glucose measuring device at least one blood glucose reference reading that is temporally associated with at least one glucose monitor data reading;  
determining a calibration equation using the at least one blood glucose reference reading  
20 and the corresponding at least one glucose monitor data reading; and  
calibrating the glucose monitor data using the calibration equation.

35. A method of calibrating characteristic monitor data, the method comprising the steps of:  
obtaining characteristic monitor data;  
obtaining from another characteristic measuring device at least one characteristic  
reference reading that is temporally associated with at least one characteristic monitor data point;  
5 calculating calibration characteristics using the at least one characteristic reference  
reading and the corresponding at least one characteristic monitor data point; and  
calibrating the obtained characteristic monitor data using the calibration characteristics.

36. A method of claim 35, wherein the at least one characteristic reference reading is at least  
10 two characteristic reference readings.

37. A method of claim 36, wherein a calculation for calculating the calibration characteristics  
is a linear regression calculation.

38. An apparatus for calibrating glucose monitor data, the apparatus comprising:  
a glucose monitor including a glucose monitor memory for storing glucose monitor data;  
a glucose sensor electronically coupled to the glucose monitor to supply the glucose  
monitor data to the glucose monitor memory;  
a blood glucose measuring device that provides at least one blood glucose reference  
20 reading that is temporally associated with at least one glucose monitor data point in the glucose  
monitor memory; and  
a processor including software to calculate calibration characteristics using the at least  
one blood glucose reference reading that is temporally associated with at least one glucose  
monitor data point and wherein the processor applies the calibration characteristics to the glucose  
25 monitor data.



39. An apparatus according to claim 38, wherein the glucose monitor includes the processor.

40. An apparatus according to claim 38, wherein the processor receives glucose monitor data from the glucose monitor.

5

41. An apparatus according to claim 38, wherein the at least one blood glucose reference reading is entered into the glucose monitor.

42. An apparatus to calibrate glucose monitor data, the apparatus comprising:

10

means for obtaining glucose monitor data;

means for obtaining from another blood glucose measuring device at least one blood glucose reference reading that is temporally associated with at least one glucose monitor data reading;

15

means for determining a calibration equation using the at least one blood glucose reference reading and the corresponding at least one glucose monitor data reading; and

means for calibrating the glucose monitor data using the calibration equation.

ABSTRACT OF THE DISCLOSURE

A method of calibrating glucose monitor data includes collecting the glucose monitor data over a period of time at predetermined intervals. It also includes obtaining at least two reference glucose values from a reference source that temporally correspond with the glucose monitor data obtained at the predetermined intervals. Also included is calculating the calibration characteristics using the reference glucose values and the corresponding glucose monitor data to regress the obtained glucose monitor data. And calibrating the obtained glucose monitor data using the calibration characteristics is included. In preferred embodiments, the reference source is a blood glucose meter, and the at least two reference glucose values are obtained from blood tests. In additional embodiments, the calculation of the calibration characteristics is obtained using linear regression and in particular embodiments, least squares linear regression. Alternatively, the calculation of the calibration characteristics is obtained using non-linear regression.



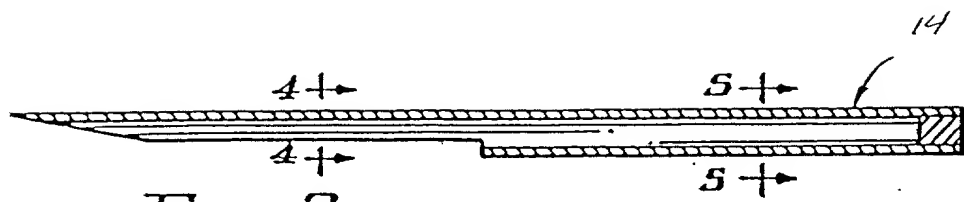


FIG. 3

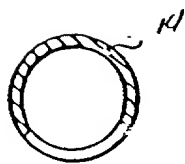


FIG. 4

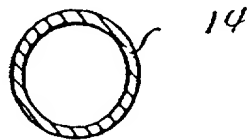


FIG. 5

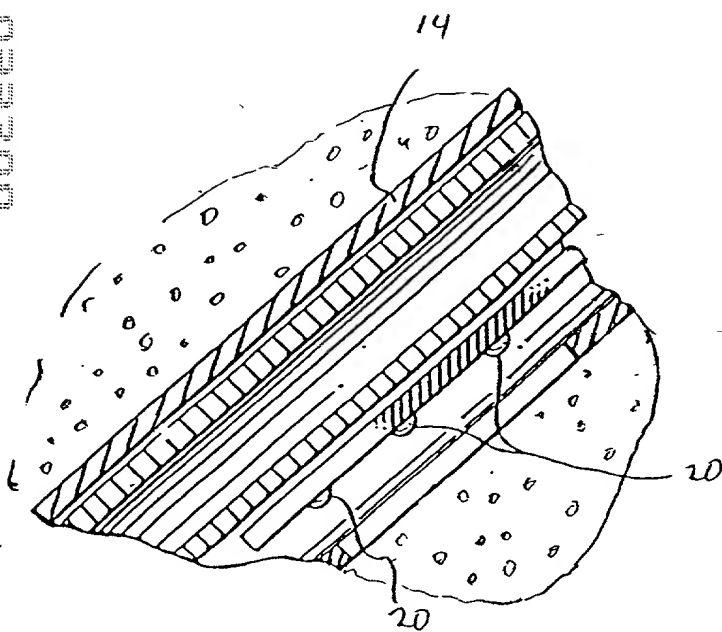


FIG. 6

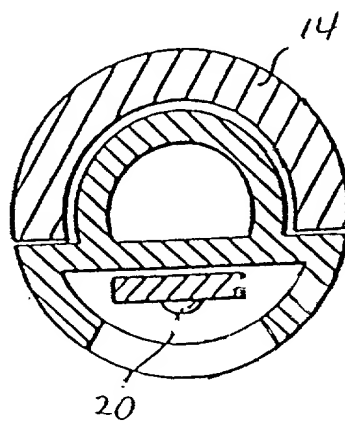


FIG. 7

2025 RELEASE UNDER E.O. 14176

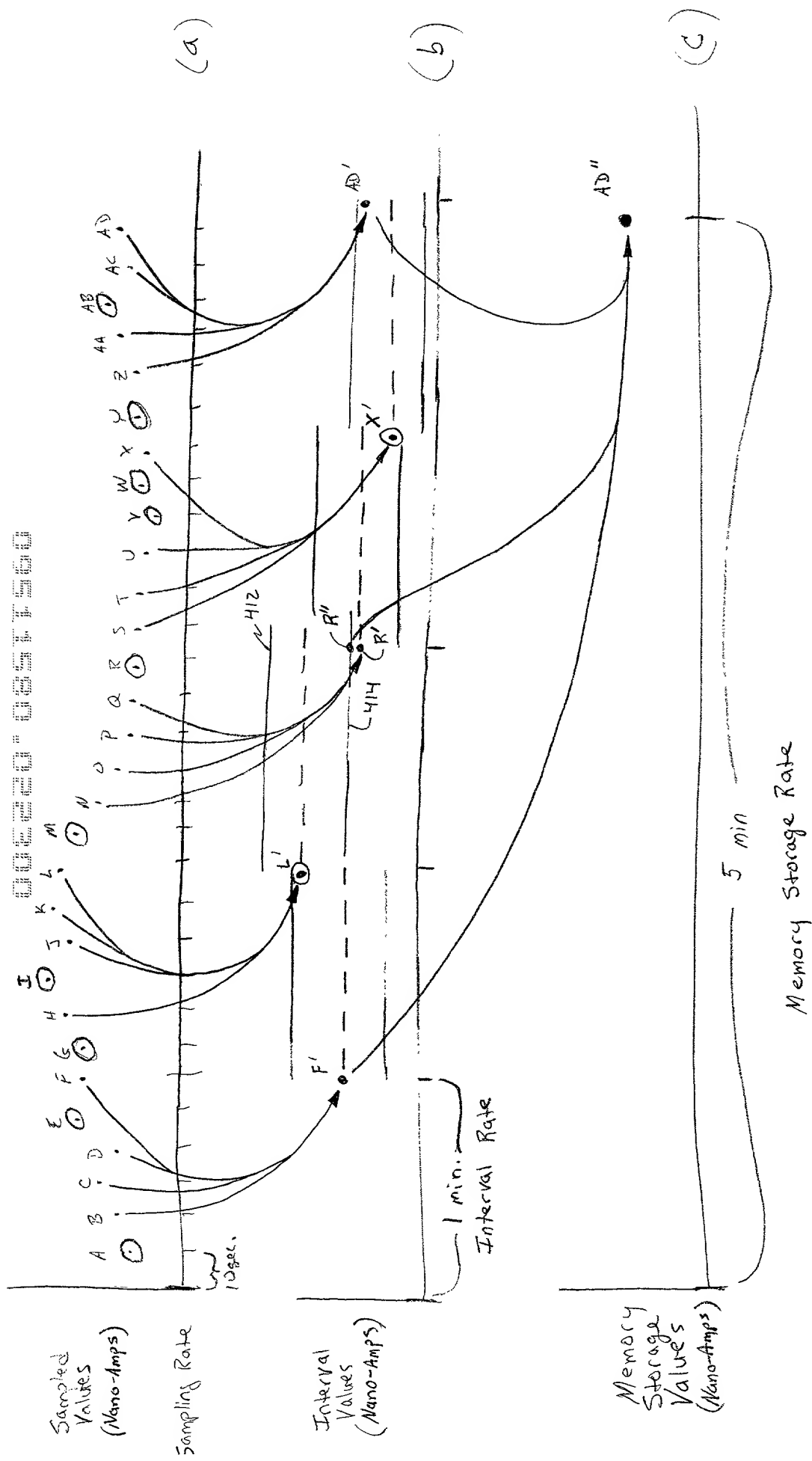


Fig. 8

NUMBER	DATE	DESCRIPTION	AMOUNT	BALANCE
1	1/1	Balance	100.00	100.00
2	1/15	Interest	1.00	101.00
3	1/31	Interest	1.00	102.00
4	2/1	Interest	1.00	103.00
5	2/15	Interest	1.00	104.00
6	2/28	Interest	1.00	105.00
7	3/1	Interest	1.00	106.00
8	3/15	Interest	1.00	107.00
9	3/31	Interest	1.00	108.00
10	4/1	Interest	1.00	109.00
11	4/15	Interest	1.00	110.00
12	4/30	Interest	1.00	111.00
13	5/1	Interest	1.00	112.00
14	5/15	Interest	1.00	113.00
15	5/31	Interest	1.00	114.00
16	6/1	Interest	1.00	115.00
17	6/15	Interest	1.00	116.00
18	6/30	Interest	1.00	117.00
19	7/1	Interest	1.00	118.00
20	7/15	Interest	1.00	119.00
21	7/31	Interest	1.00	120.00
22	8/1	Interest	1.00	121.00
23	8/15	Interest	1.00	122.00
24	8/31	Interest	1.00	123.00
25	9/1	Interest	1.00	124.00
26	9/15	Interest	1.00	125.00
27	9/30	Interest	1.00	126.00
28	10/1	Interest	1.00	127.00
29	10/15	Interest	1.00	128.00
30	10/31	Interest	1.00	129.00
31	11/1	Interest	1.00	130.00
32	11/15	Interest	1.00	131.00
33	11/30	Interest	1.00	132.00
34	12/1	Interest	1.00	133.00
35	12/15	Interest	1.00	134.00
36	12/31	Interest	1.00	135.00
37	1/1	Interest	1.00	136.00
38	1/15	Interest	1.00	137.00
39	1/31	Interest	1.00	138.00
40	2/1	Interest	1.00	139.00
41	2/15	Interest	1.00	140.00
42	2/28	Interest	1.00	141.00
43	3/1	Interest	1.00	142.00
44	3/15	Interest	1.00	143.00
45	3/31	Interest	1.00	144.00
46	4/1	Interest	1.00	145.00
47	4/15	Interest	1.00	146.00
48	4/30	Interest	1.00	147.00
49	5/1	Interest	1.00	148.00
50	5/15	Interest	1.00	149.00
51	5/31	Interest	1.00	150.00
52	6/1	Interest	1.00	151.00
53	6/15	Interest	1.00	152.00
54	6/30	Interest	1.00	153.00
55	7/1	Interest	1.00	154.00
56	7/15	Interest	1.00	155.00
57	7/31	Interest	1.00	156.00
58	8/1	Interest	1.00	157.00
59	8/15	Interest	1.00	158.00
60	8/31	Interest	1.00	159.00
61	9/1	Interest	1.00	160.00
62	9/15	Interest	1.00	161.00
63	9/30	Interest	1.00	162.00
64	10/1	Interest	1.00	163.00
65	10/15	Interest	1.00	164.00
66	10/31	Interest	1.00	165.00
67	11/1	Interest	1.00	166.00
68	11/15	Interest	1.00	167.00
69	11/30	Interest	1.00	168.00
70	12/1	Interest	1.00	169.00
71	12/15	Interest	1.00	170.00
72	12/31	Interest	1.00	171.00
73	1/1	Interest	1.00	172.00
74	1/15	Interest	1.00	173.00
75	1/31	Interest	1.00	174.00

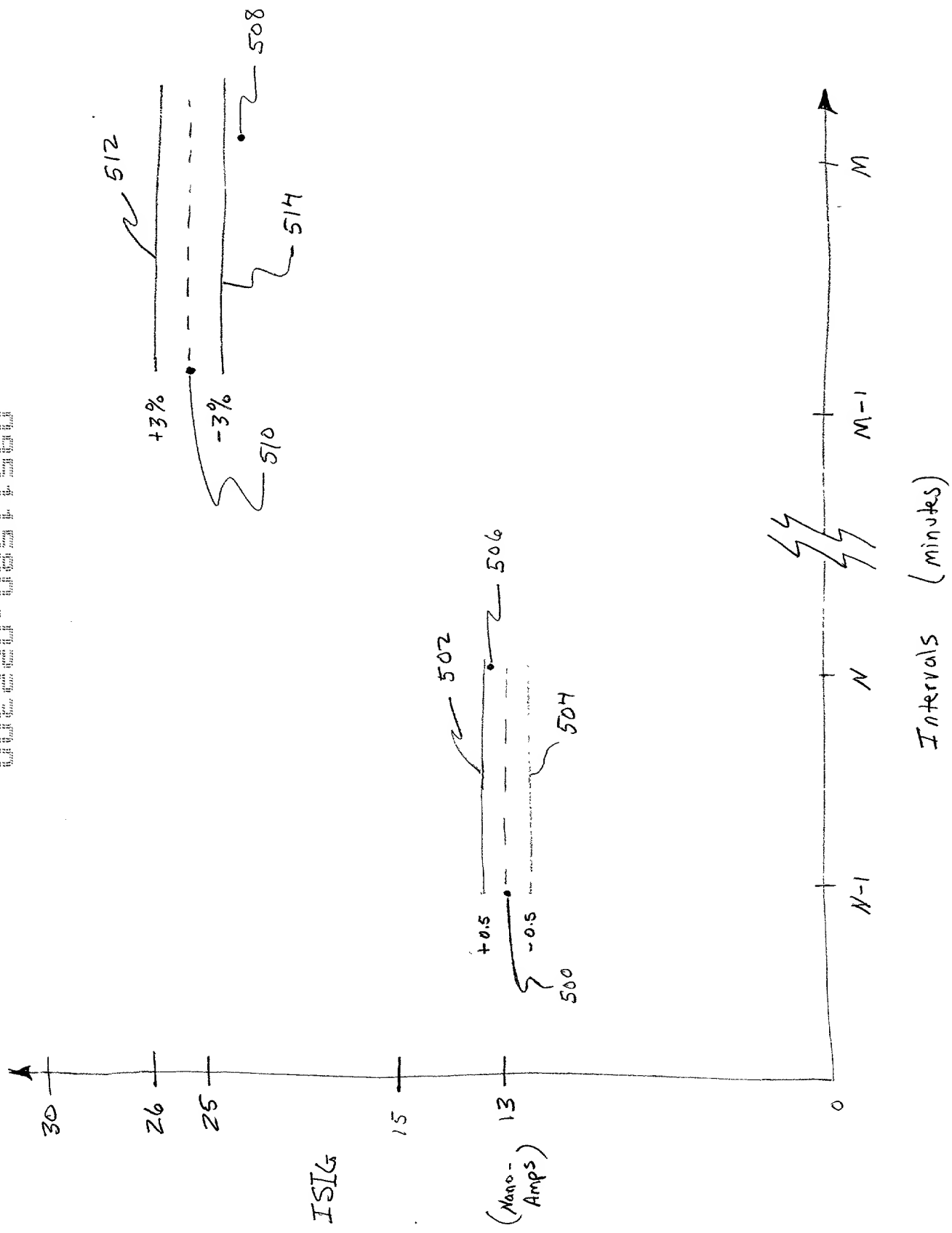


Fig. 9



Species	1990-1991	1991-1992	1992-1993	1993-1994	1994-1995	1995-1996	1996-1997	1997-1998	1998-1999	1999-2000	2000-2001	2001-2002	2002-2003	2003-2004	2004-2005	2005-2006	2006-2007	2007-2008	2008-2009	2009-2010	2010-2011	2011-2012	2012-2013	2013-2014	2014-2015	2015-2016	2016-2017	2017-2018	2018-2019	2019-2020	2020-2021	2021-2022	2022-2023	2023-2024	2024-2025	2025-2026	2026-2027	2027-2028	2028-2029	2029-2030	2030-2031	2031-2032	2032-2033	2033-2034	2034-2035	2035-2036	2036-2037	2037-2038	2038-2039	2039-2040	2040-2041	2041-2042	2042-2043	2043-2044	2044-2045	2045-2046	2046-2047	2047-2048	2048-2049	2049-2050	2050-2051	2051-2052	2052-2053	2053-2054	2054-2055	2055-2056	2056-2057	2057-2058	2058-2059	2059-2060	2060-2061	2061-2062	2062-2063	2063-2064	2064-2065	2065-2066	2066-2067	2067-2068	2068-2069	2069-2070	2070-2071	2071-2072	2072-2073	2073-2074	2074-2075	2075-2076	2076-2077	2077-2078	2078-2079	2079-2080	2080-2081	2081-2082	2082-2083	2083-2084	2084-2085	2085-2086	2086-2087	2087-2088	2088-2089	2089-2090	2090-2091	2091-2092	2092-2093	2093-2094	2094-2095	2095-2096	2096-2097	2097-2098	2098-2099	2099-2100	2100-2101	2101-2102	2102-2103	2103-2104	2104-2105	2105-2106	2106-2107	2107-2108	2108-2109	2109-2110	2110-2111	2111-2112	2112-2113	2113-2114	2114-2115	2115-2116	2116-2117	2117-2118	2118-2119	2119-2120	2120-2121	2121-2122	2122-2123	2123-2124	2124-2125	2125-2126	2126-2127	2127-2128	2128-2129	2129-2130	2130-2131	2131-2132	2132-2133	2133-2134	2134-2135	2135-2136	2136-2137	2137-2138	2138-2139	2139-2140	2140-2141	2141-2142	2142-2143	2143-2144	2144-2145	2145-2146	2146-2147	2147-2148	2148-2149	2149-2150	2150-2151	2151-2152	2152-2153	2153-2154	2154-2155	2155-2156	2156-2157	2157-2158	2158-2159	2159-2160	2160-2161	2161-2162	2162-2163	2163-2164	2164-2165	2165-2166	2166-2167	2167-2168	2168-2169	2169-2170	2170-2171	2171-2172	2172-2173	2173-2174	2174-2175	2175-2176	2176-2177	2177-2178	2178-2179	2179-2180	2180-2181	2181-2182	2182-2183	2183-2184	2184-2185	2185-2186	2186-2187	2187-2188	2188-2189	2189-2190	2190-2191	2191-2192	2192-2193	2193-2194	2194-2195	2195-2196	2196-2197	2197-2198	2198-2199	2199-2200	2200-2201	2201-2202	2202-2203	2203-2204	2204-2205	2205-2206	2206-2207	2207-2208	2208-2209	2209-2210	2210-2211	2211-2212	2212-2213	2213-2214	2214-2215	2215-2216	2216-2217	2217-2218	2218-2219	2219-2220	2220-2221	2221-2222	2222-2223	2223-2224	2224-2225	2225-2226	2226-2227	2227-2228	2228-2229	2229-2230	2230-2231	2231-2232	2232-2233	2233-2234	2234-2235	2235-2236	2236-2237	2237-2238	2238-2239	2239-2240	2240-2241	2241-2242	2242-2243	2243-2244	2244-2245	2245-2246	2246-2247	2247-2248	2248-2249	2249-2250	2250-2251	2251-2252	2252-2253	2253-2254	2254-2255	2255-2256	2256-2257	2257-2258	2258-2259	2259-2260	2260-2261	2261-2262</
---------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-----------	-------------

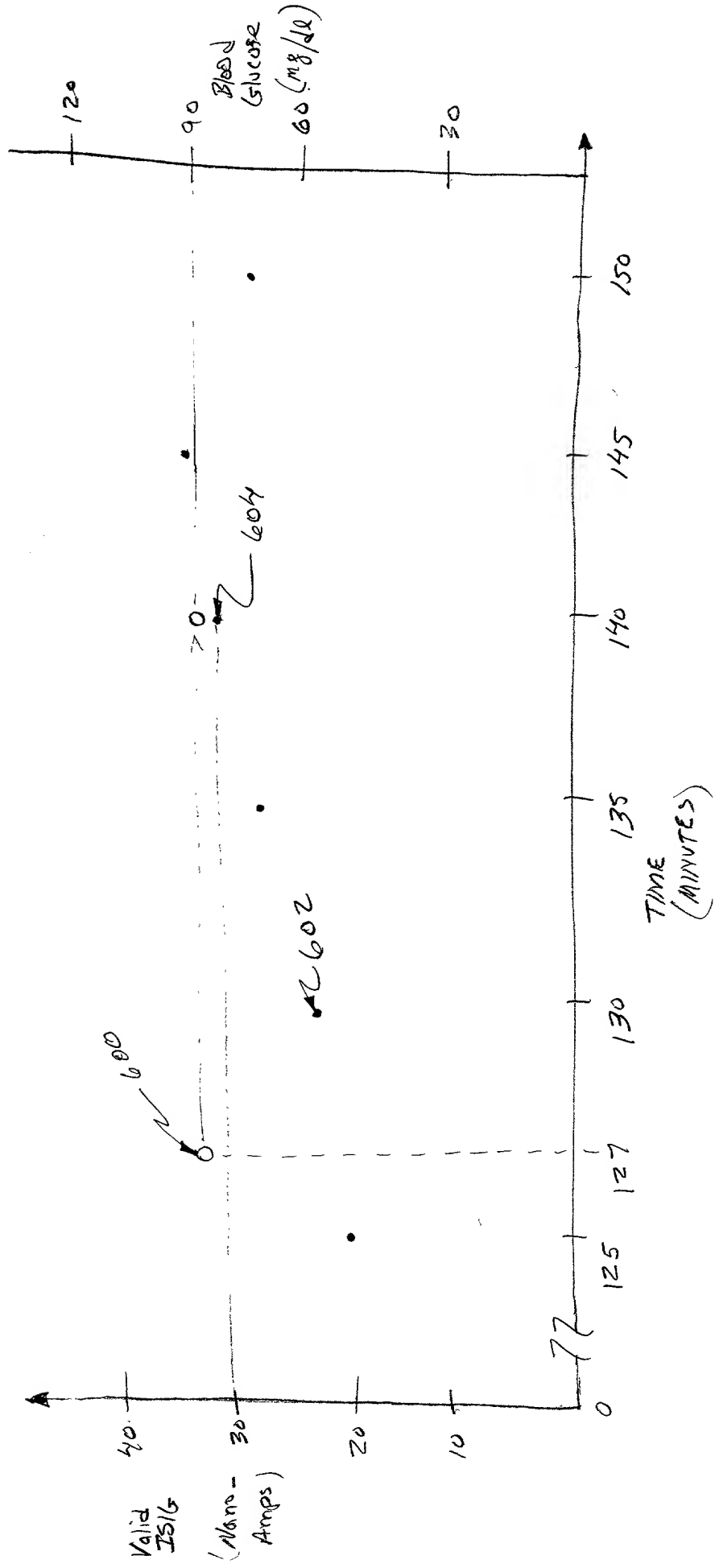


Fig. 11



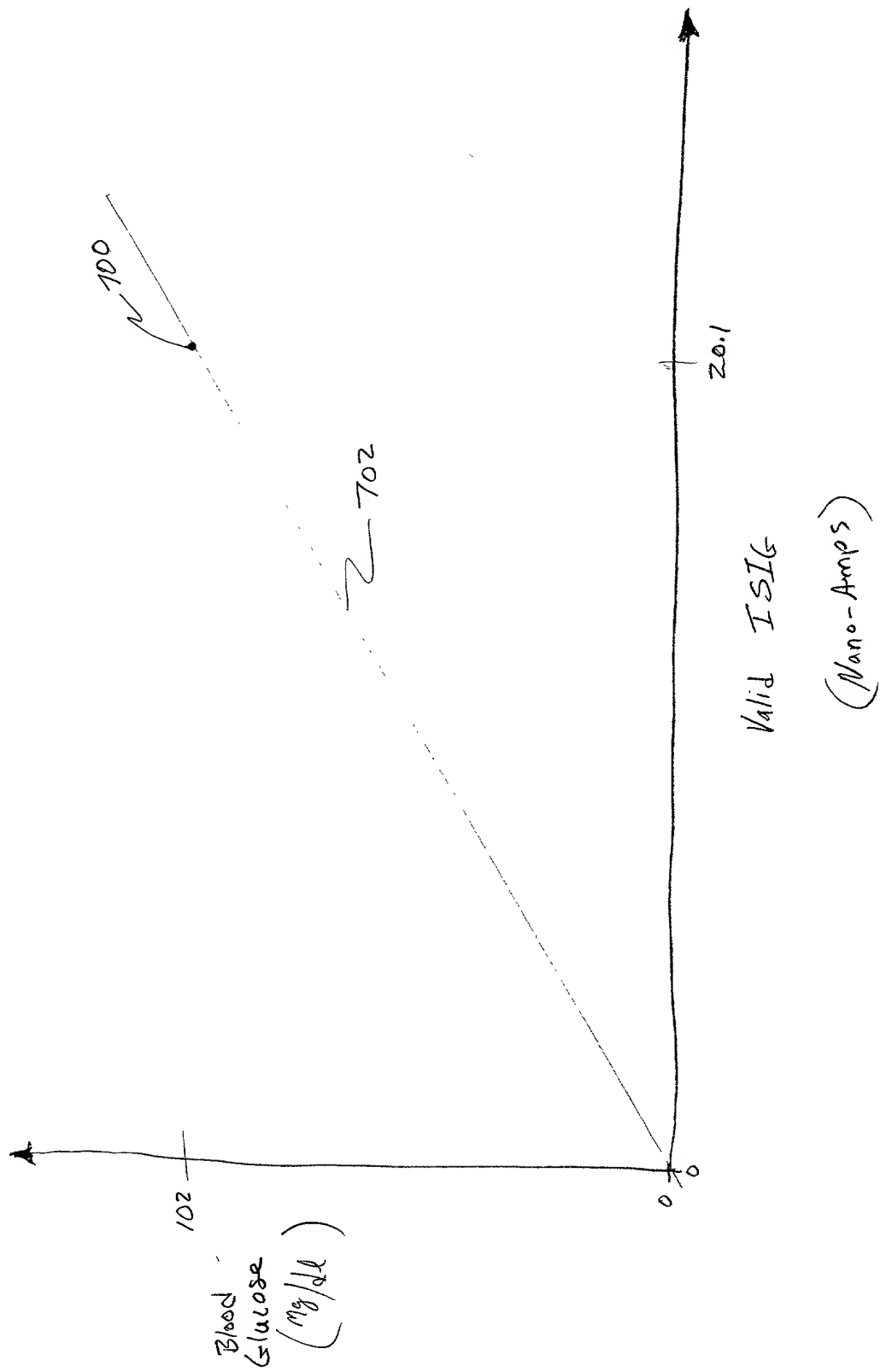


Fig. 12

Initialization / Standardization  
Complete

Pair a Blood Glucose Standard Reading  
with a Valid ISIG Memory Storage Value.

Calculate Single Point Sensitivity Ratio (SPSR)  
$$\text{SPSR} = \frac{\text{Blood Glucose Standard Reading}}{\text{Valid ISIG}}$$

Select Offset Value  
 $\text{SPSR} < 7 \Rightarrow \text{offset} = 3$   
 $\text{SPSR} \geq 7 \Rightarrow \text{offset} = 0$

Calculate Modified SPSR (MSPSR)  
$$\text{MSPSR} = \frac{\text{Blood Glucose Standard Reading}}{(\text{Valid ISIG} - \text{offset})}$$

Initial Calibration is complete

$$\text{Blood Glucose} = (\text{Valid ISIG} - \text{offset}) \times \text{MSPSR}$$

Fig. 13

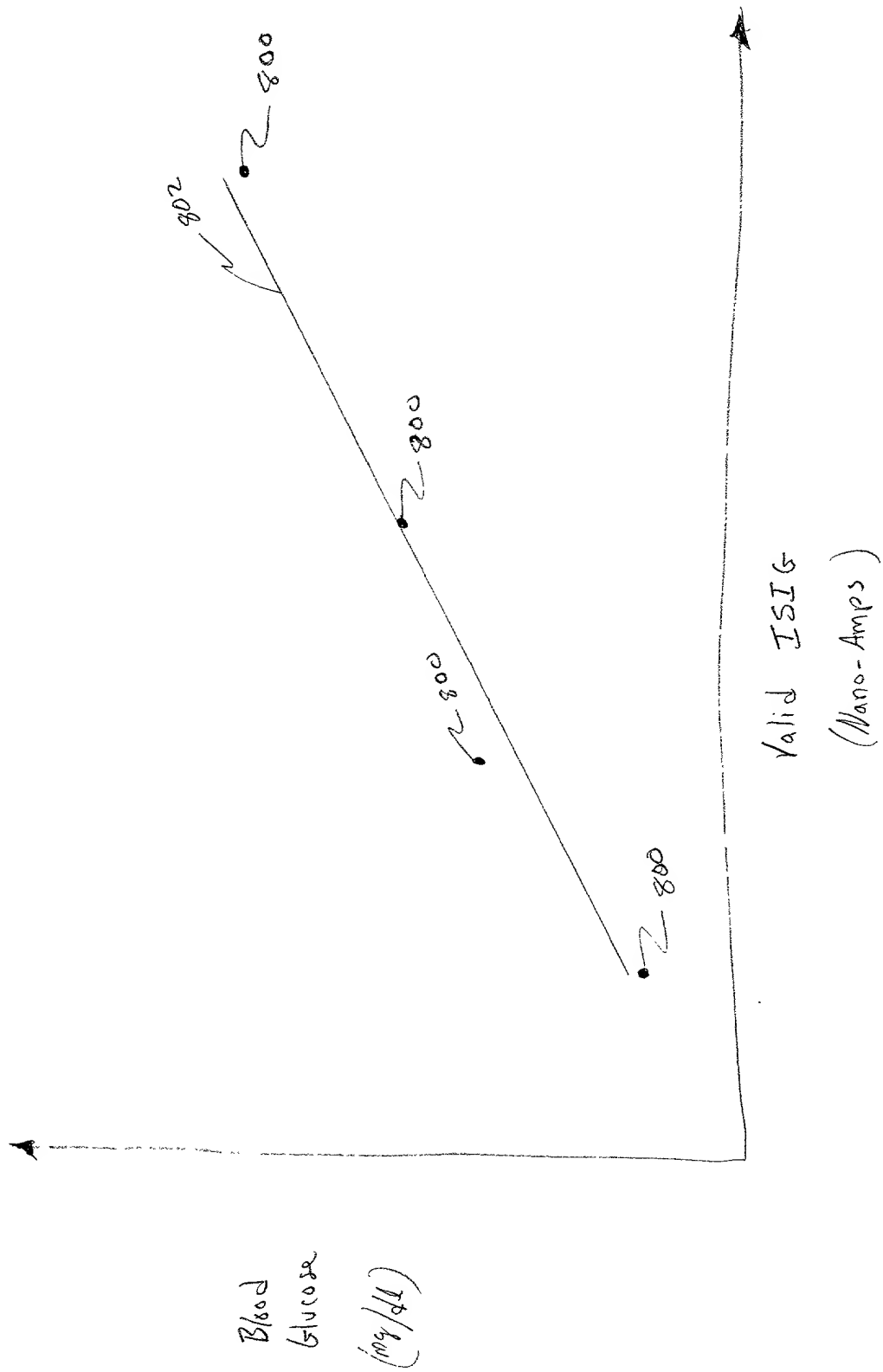
[illegible]

Fig. 14

2 or more paired calibration data points are available.

Calculate Linear Regression Sensitivity Ratio (LRSR)

$$LRSR = \frac{\sum_{i=1}^N [X_i Y_i]}{\sum_{i=1}^N [X_i^2]}$$

$i$  = the  $i$ th pair  
where  $X$  = Valid ISIG  
 $Y$  = Blood Glucose  
Standard  
Reading

$N$  = Number of paired  
calibration data points

Select Offset Value

$LRSR < 7 \Rightarrow \text{offset} = 3$

$LRSR \geq 7 \Rightarrow \text{offset} = 0$

Calculate Modified LRSR (MLRSR)

$$MLRSR = \frac{\sum_{i=1}^N [(X_i - \text{offset}) Y_i]}{\sum_{i=1}^N [(X_i - \text{offset})^2]}$$

Calibration is complete

Blood Glucose = (Valid ISIG - offset) (MLRSR)

Fig. 15